

Large Volume Centrifuges

■ Z496-K Floor-Standing

■ Z496-K-UC Underbench



HERMLE
Benchmark



■ The Perfect Fit



Easy Handling:

- The ergonomic control panel is very user-friendly
- No need to remove gloves during operation
- Strictly a HERMLE application, you can store up to 99 different programs.
- Both versions are equipped with wheels, to easily maneuver around the laboratory, to the most convenient location. There is also a mechanism to place them on the feet as well
- Make yourself comfortable; operate this model from your laptop or computer. Our multi-functional software, "H-Control", allows you to control the centrifuge independently and save your protocols, (this function has to be ordered separately).
- An optional alarm will sound when your run is finished

Safety:

- The Z496-K has an automatic imbalance shut off system.
- There is over speed protection incorporated for each rotor.
- The centrifuge will not start until the lid is locked.
- Error code system to quickly detect issues
- Stable and firm design, made of steel, ensures excellent mechanical reliability and protection.

Efficient and Versatile:

- Due to the wide range of bucket and rotor options, the Z496 can hold all common tubes and bottles, up to 1,000 ml.
- "Pre-cooling" function specific to each rotor for maximum efficiency
- Keep dangerous and hazardous aerosols locked; hermetically sealed buckets available
- Save space in your laboratory. The Underbench model is easily placed below a standard laboratory table.
- The Floor-Standing mode; is supplied with a storage drawer which is very useful for your accessories.

4 x 1,000 ml

Z496-K



**Underbench
Centrifuge**



**Floor-Standing
Centrifuge**

Features:

- Maintenance-free induction drive
- Microprocessor with large LCD display
- Motor driven lid lock
- Automatic radio-frequency rotor identification system (RFID) with over speed protection, immediately upon rotor insertion
- Active imbalance identification and cut-off (rotor-specific)
- CFC-free refrigeration system
- Wide range of rotors and accessories
- Audible signal at the end of each run, with several melody options
- Stainless steel chamber
- Manufactured in accordance to International Safety Regulations, i.e. IEC 61010

The Z496-K centrifuges offer the highest tube/bottle capacity in our range. Several rotors are available, including a popular 4 position swinging bucket rotor with a wide offering of accessories for the centrifugation of tubes, bottles, plates and blood bags. The high capacity rectangular buckets accept microplates directly or a variety of tube inserts for accommodating as many as 40 x 50ml conical or 96 x 15ml conical tubes per run. A standard floor standing model is available (height: 38.6 in.) which includes a convenient lower drawer for the storage of frequently used rotors and accessories. In addition, a space saving "underbench" model (height: 27.5 in.) is available when under the counter storage is preferred. Both models are equipped with wheels for easy maneuverability around the laboratory. For all your high capacity, high volume needs, the Z496-K floor standing series is the perfect centrifuge.

Technical Data

Max. Speed:	10,500 rpm
Max. RCF:	16,022 xg
Max. Volume:	4 x 1,000 ml
Speed range:	200 - 10,500 rpm
Running time:	59 min 50 sec / 10 s increments 99 h 59 min / 1 min increments

Dimensions (W x H x D):

Underbench Unit:	24.4 in x 27.5 in x 27 in 62 cm x 70 cm x 69 cm
Floor-Standing Unit:	24.4 in x 38.5 in x 27 in 62 cm x 98 cm x 69 cm

Weight:

Underbench Unit:	346 lbs / 157 kg
Floor-Standing Unit:	419 lbs / 190 kg

Power input: 2.3 KW

Ordering Information

Underbench: Z496-K-UC, 230 V / 50 - 60Hz

Floor-Standing: Z496-K, 230 V / 50 - 60 Hz

Interface incl. operating software: Z496-HCTRL



Distinct Control Panel

- Simple one-handed operation
- Easy to program with gloves on
- Splash-proof
- Permanent indication of pre-set and actual values
- Selection of speed in both rpm and g-force, in increments of 10
- 10 acceleration and deceleration rates, possibility of unbraked deceleration
- Temperature range from -20°C to 40°C, in increments of 1°C
- Pre-cooling program
- Standstill cooling
- Pre-selection of the running time from 10 s to 99 h 59 min or continuous
- Storage of up to 99 runs, including rotor
- Quick-key for short runs
- Timer: the centrifuge starts at the preselected time

Accessories Z496-K

Rotor/Bucket

Order No.:
Z496-1000-AC*
4 x 1000ml Swing out rotor
Angle: 0-90°



*Z496-1000
Swing out rotor with
rectangular buckets included

Adapter

Order No.:
Z496-1000-1L
Round bucket



Acceleration in s: 84/328
Deceleration in s: 50/1,067
Max. Speed: 4,600 rpm
Max. RCF: 5,203 xg
Max. Radius: 22 cm
Temp. at Max. Speed: 9°C



Z446-750-LID
PC-Screw cap hermetical
sealable, for round bucket
Z496-1000-1L

Order No.:
Z496-1000-RBUK
Rectangular bucket/Microplate bucket










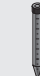





Acceleration in s: 84/302
Deceleration in s: 53/1,178
Max. Speed: 4,200 rpm
Max. RCF: 3,964 xg
Max. Radius: 20.1 cm
Max. Loading Height: 92 mm
Max. Loading Weight: 1,200 g
Temp. at Max. Speed: 4°C



Z496-1000-MPRA
Removing aid for
microplates

Adapters for Rectangular Bucket

		Corning				Falcon						
Rectangular Bucket Order No.: Z496-1000-RBUK	Volume	750 ml	650 ml	500 ml*	500 ml	250/290 ml	175/225**ml	100 ml	50 ml	50 ml	15 ml	5 ml
	ø in mm	99	99	97	99	62	61	40	29	29	17	17
												
	Length min./max. (mm)	152	136	148	114	122/137	118/137	102/123	90/120	116	120	60/70
	Tubes per Rotor	4	4	4	4	8	8	12	40	40	96	96
	Radius Correction (cm)	0.5	0.5	0.5	0.5	0.4	0.5	0.5	0.1	0.1	0.1	0.1
	Order No. Z496-1000_	-A750	-A750	-A500	-A500	-A250	-A225	-A100	-A50R	-A50	-A15	-A20

Vacutainer®/Sarstedt®		
Volume	5.5-12 ml	2.6-8 ml
ø in mm	16	13
		
Length min./max. (mm)	60/120	65/100
Tubes per Rotor	152	192
Radius Correction (cm)	0.1	0.1
Order No. Adapter	-A10	-A57

*Only applicable with combination of Z496-1000-A750 + Z446-750-A500

Adapters for Round Bucket

Volume	Corning®					Falcon®						
	1000 ml	750 ml	650 ml	500 ml	500 ml	250 ml	250/290 ml	175/225 ml	50 ml	50 ml	30 ml	
Round Bucket Order No.: Z496-1000-1L	99	99	99	97	99	69.5	60	62	61	29	29	25
ø in mm												
Length min./max. (mm)	173/178	152	136	148	114	165	159	122/137	118/137	102/123	116	102/123
Tubes per Rotor	4	4	4	4	4	4	4	4	4	28	20	32
Radius Correction (cm)	0	0	0	0.5	0	0.5	0.4	0.5	0.2	0.1	0.1	0.1
Order No. Z496-1L_	without	without	without	*-A500	without	-A500S	-A250F	*-A250	*-A225	*-A50R	*-A50-E	*-A30

Vacutainer®			
Volume	12/15 ml	15 ml	6-10 ml
ø in mm	17	17	16
Length min./max. (mm)	97/120	120	75/100
Tubes per Rotor	72	48	72
Radius Correction (cm)	0.1	0.1	0.3
Order No. Adapter	*-A15R	-A15	*-A10

*use Z446-750 as prefix

Accessories Z496-K

Rotor/Bucket

Order No.: Z496-1000-AC

4 x 1000ml Swing out rotor
Angle: 0-90°



Adapter

Order No.: Z496-500-BLD

Blood bag bucket for 500 ml



Acceleration in s: 62/328
Deceleration in s: 48/981
Max. Speed: 4,600 rpm
Max. RCF: 5,298 xg
Max. Radius: 22.4 cm
Temp. at Max. Speed: 7°C



Z496-500-ADP

Adapter for bucket
Z496-500-ADP

Order No.: Z496-750-BLD

Blood bag bucket for 750 ml



Acceleration in s: 82/324
Deceleration in s: 50/1,111
Max. Speed: 4,600 rpm
Max. RCF: 5,274 xg
Max. Radius: 22.3 cm
Temp. at Max. Speed: 9°C



Z496-750-ADP

Adapter for bucket
Z496-750-ADP

Order No.:
Z446-04500



Angle rotor 4 x 500 ml
Angle: 24°

Acceleration in s: 73/575
Deceleration in s: 84/2,317
max. speed: 8,000 rpm
max. RCF: 10,367 xg
max. radius: 14.5 cm
max. ø: 69.5 mm
Temp. at max. speed: 4°C

Order No.:

without Adapter
Z446-500-A250F
Z446-500-A250
Z446-500-A225
Z446-500-A225
Z446-500-A50
Z446-500-A15C

Tube rack

1 x 500 ml
1 x 250 ml
1 x 250 ml conical
1 x 175 ml Falcon®
1 x 225 ml Falcon®
1 x 50 ml conical
6 x 15 ml conical

ø (mm) Length (min./max.)

69.5 165
62 122
60 159
61 118
61 137
29 116
17 120

Order No.:
Z366-06250



Angle rotor 6 x 250 ml
Angle: 28°

Acceleration in s: 66/573
Deceleration in s: 84/2,317
max. speed: 8,000 rpm
max. RCF: 10,016 xg
max. radius: 14 cm
max. ø: 62 mm
Temp. at max. speed: -1°C

Order No.:

without Adapter
Z366-06250-A50
Z366-06250-A50R
Z366-06250-A30
Z366-06250-A15

Tube rack

1 x 250 ml
1 x 50 ml conical
1 x 50 ml
2 x 30 ml
4 x 15 ml conical

Ø (mm) Length (min./max.)

62 122
29 118
29 103
26 92
17 120

Order No.:
Z326-1050



Angle rotor 10 x 50 ml conical
Angle: 38°

Acceleration in s: 53/441
Deceleration in s: 71/1,411
Max. Speed: 10,500 rpm
Max. RCF: 16,022 xg
Max. Radius: 13.0 cm
Max. ø: 29 mm
Temp. at max. speed: -1°C

Order No.:

Adapter
Z326-50-B15HS 1 x 15 ml RB
Z326-15-A7* 1 x 7 ml
Z326-15-A5* 1 x 5 ml
Z326-50-A15 1 x 15 ml conical

ø (mm) Length (min./max.)

17 97/123
13.5 95/115
13.5 70/90
17 120

*only possible with adapter Z326-50-A15 or Z326-50-B15

Z496-K Floor-Standing Model

The Floor-Standing Model, features a storage drawer for easy access to your accessories.



www.HermleUSA.com

PH: 908-769-5555

EM: info@benchmarkscientific.com

Distributed by:

NELSON JAMESON
INC.

800-826-8302 nelsonjameson.com



Technical modification rights reserved. V 10/17